



Mini-Monobloc Fill/Finish Packaging System

General Overview

- Ideal for small batch production and utilizing one operator
- Small footprint includes sorting bowls and controls enclosure
- Designed for HEPA filtration unit integration
- Easy mobility with lockable castors
- Configurable for fill/cap, fill/tip/cap, & fill/stopper
- User-friendly, menu driven HMI touch screen operator interface
- Runs plastic & glass bottles, vials, microtubes, and pre-capped microtubes



Unique Performance. Undeniable Quality. Unrivaled Support.
Go Further, Do More with Filamatic.

***Mini-Monobloc Fill/Finish
Packaging System***


866.258.1914
info@filamatic.com
www.filamatic.com




Technical Features

- Approximate dimensions: 24" x 48" x 54"
- Fill accuracy: + / - 0.5% or better
- Fill volume range: microliters to 120ml
- Choice of filling system: positive displacement, peristaltic, or cicro-dose rotary piston
- Designed with casters for portability


Handfeed - Fill - Tip - Screw Cap

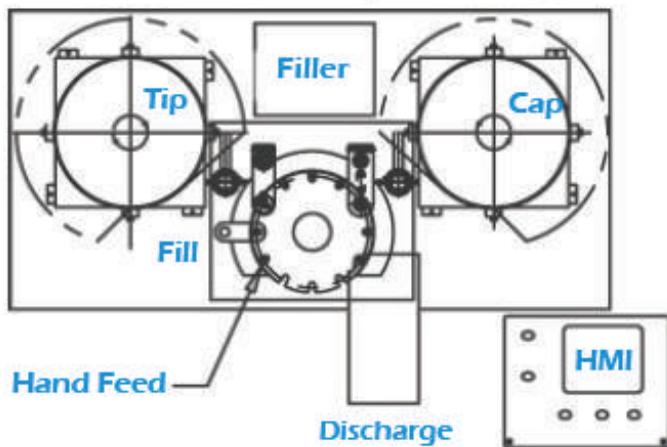


Feed - Fill - Cap (Microtubes)

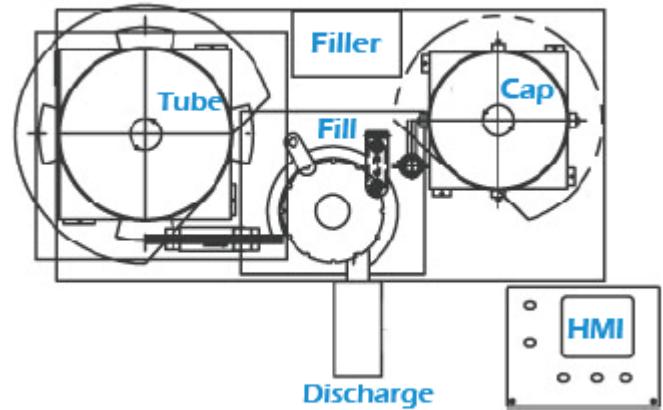

Fill - Stopper - Crimp - (Vial)

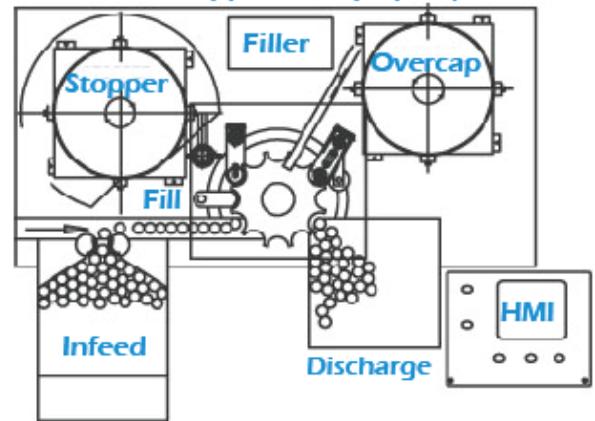

Feed - Uncap - Fill - Recap
(Precapped Microtubes)

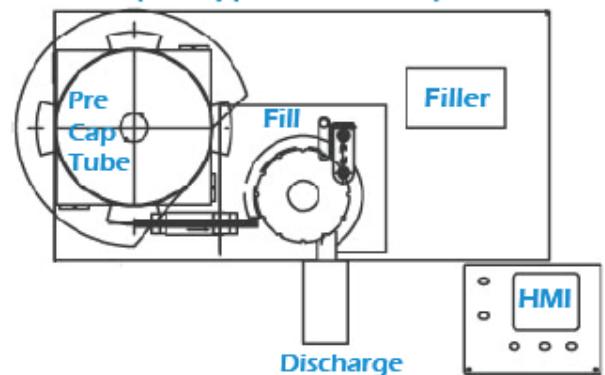

Filamatic

4119 Fordleigh Road • Baltimore, MD 21215
866.258.1914 • info@filamatic.com • www.filamatic.com